



# SURGICRYL<sup>®</sup>

910


POLYGLACTINE

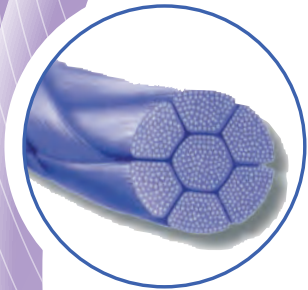

<b>Type</b>	Braided and coated multifilament
<b>Composition</b>	Polyglactine 910, a copolymer made of 90% glycolide and 10%L-lactide
<b>Coating</b>	Poly-glycolide-co-L-lactide and calcium stearate (<1%)
<b>Colour</b>	Violet
<b>Tissue reaction</b>	Minimal
<b>Absorption</b>	The hydrolytic action by which the material is broken down results in total absorption between 56 and 70 days. Between 40% and 50% of tensile strength remain after 21 days.
<b>Presentation</b>	Box with 12 sutures

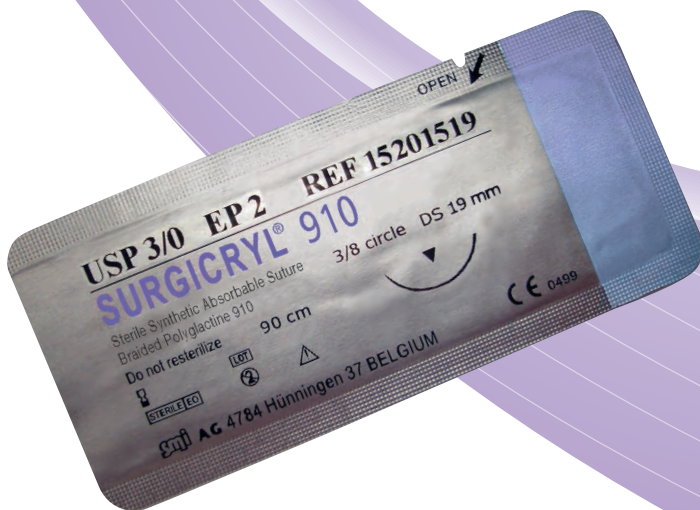

*For quality  
&  
safety*

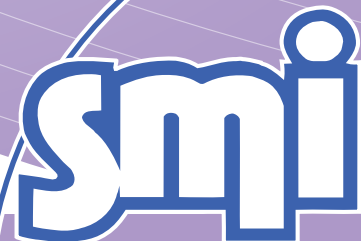

[www.sutures.be](http://www.sutures.be)

# SURGICRYL®

## 910

### POLYGLACTINE

### CHARACTERISTICS

- High tensile strength
- Synthetic
- Breaks down by hydrolysis
- Braided structure
  
- Coated multifilament
  
- Hermetically sealed packing

### BENEFITS

- Excellent tissue support
- Minimal tissue reaction
- Predictable absorption rate
- Excellent handling properties
  - High flexibility
  - High tensile strength
  - Excellent knot security
- Soft passage through the tissues with minimal sawing, tissue drag and trauma
- Good knot tie down / adjustability
- Reduced capillary action
  
- Guaranteed seal and product sterility

**INDICATIONS** General soft tissue closing and/or ligation especially in


- General surgery
- Skin closure
- Gastrointestinal surgery
- Plastic surgery
- Ophthalmic surgery
- Gynaecology
- Obstetrics
- Orthopaedics
- Urology

**CONTRAINDICATIONS** Cardiovascular surgery  
Neurosurgery


Being absorbable not to use when an extended wound support is required.



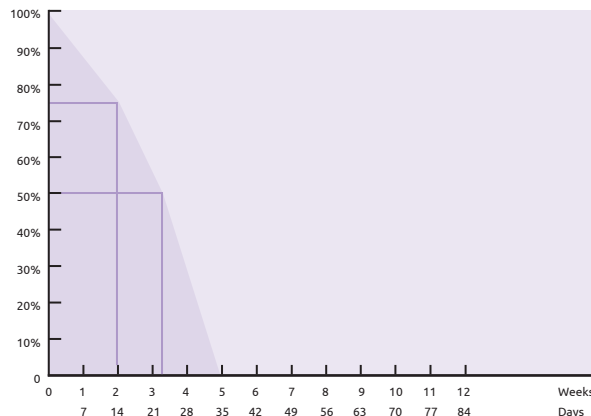


Caution in patients with delayed wound healing (e.g. compromised immune system).



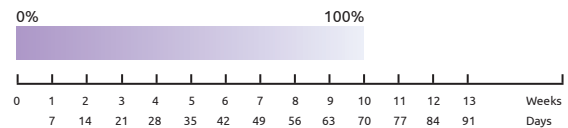

### TENSILE STRENGTH




















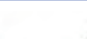
### ABSORPTION


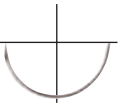








Round bodied taper point			Metric	USP	Length	Order number
1/2 Circle 	HR 13 mm		0.7	6/0	75 cm	<i>violet</i> 15070113
			1	5/0	75 cm	<i>violet</i> 15100113
			1.5	4/0	75 cm	<i>violet</i> 15150113
	HR 17 mm		1	5/0	75 cm	<i>violet</i> 15100117
			1.5	4/0	75 cm	<i>violet</i> 15150117
			2	3/0	75 cm	<i>violet</i> 15200117
	HR 22 mm		3	2/0	75 cm	<i>violet</i> 15300117
			1.5	4/0	75 cm	<i>violet</i> 15150122
			2	3/0	75 cm	<i>violet</i> 15200122
			3	2/0	75 cm	<i>violet</i> 15300122

Round bodied taper point			Metric	USP	Length	Order number
3/8 Circle 	DR 17 mm		1.5	4/0	45 cm	<i>violet</i> 15150217
	DR 22 mm		2	3/0	75 cm	<i>violet</i> 15200222

Reverse cutting			Metric	USP	Length	Order number
1/2 Circle 	HS 12 mm		0.7	6/0	45 cm	<i>violet</i> 15071412
	HS 17 mm		1	5/0	75 cm	<i>violet</i> 15101417
			1.5	4/0	75 cm	<i>violet</i> 15151417
	HS 22 mm		2	3/0	75 cm	<i>violet</i> 15201422
	HS 23 mm		5	2	90 cm	<i>violet</i> 15501423
	HS 26 mm		2	3/0	75 cm	<i>violet</i> 15201426
3			2/0	75 cm	<i>violet</i> 15301426	

Reverse cutting			Metric	USP	Length	Order number	
3/8 Circle 	DS 12 mm		0.7	6/0	45 cm	<i>violet</i> 15071512	
			1	5/0	45 cm	<i>violet</i> 15101512	
			1.5	4/0	45 cm	<i>violet</i> 15151512	
	DS 16 mm		0.7	6/0	45 cm	<i>violet</i> 15071516	
			1	5/0	75 cm	<i>violet</i> 15101516	
			1.5	4/0	75 cm	<i>violet</i> 15151516	
				2	3/0	45 cm	<i>violet</i> 15201516-45
	DS 16 mm Easy Pass		2	3/0	50 cm	<i>undyed</i> 15201516UP	
	DS 16 mm black		0.7	6/0	45 cm	<i>violet</i> 15071516VB	
	DS 19 mm		1	5/0	75 cm	<i>violet</i> 15101519	
			1.5	4/0	75 cm	<i>violet</i> 15151519	
			2	3/0	75 cm	<i>violet</i> 15201519	
DS 19 mm Easy Pass		1.5	4/0	50 cm	<i>undyed</i> 15151519UP		
DS 19 mm black		1	5/0	75 cm	<i>violet</i> 15101519VB		
		1.5	4/0	75 cm	<i>violet</i> 15151519VB		
DS 24 mm		1.5	4/0	75 cm	<i>violet</i> 15151524		
		2	3/0	75 cm	<i>violet</i> 15201524		
		3	2/0	75 cm	<i>violet</i> 15301524		
DS 24 mm Easy Pass		2	3/0	50 cm	<i>undyed</i> 15201524UP		

Round bodied taper cutting point			Metric	USP	Length	Order number
1/2 Circle 	HRT 13 mm		1.5	4/0	75 cm	<i>violet</i> 15150513
	HRT 17 mm		1.5	4/0	75 cm	<i>violet</i> 15150517
			2	3/0	75 cm	<i>violet</i> 15200517

Round bodied taper cutting point			Metric	USP	Length	Order number
3/8 Circle 	DRT 13 mm		1.5	4/0	75 cm	<i>violet</i> 15150413
	DRT 17 mm		1.5	4/0	75 cm	<i>violet</i> 15150417
			2	3/0	75 cm	<i>violet</i> 15200417