


SURGICRYL[®]

MONOFILAMENT

POLYDIOXANONE

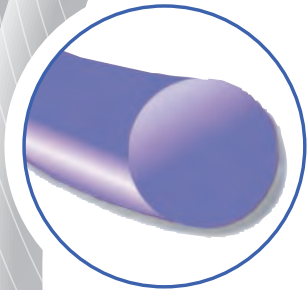

Type	Monofilament
Composition	Polydioxanone, a polymer made from polyester poly (p-dioxanone)
Coating	None
Colour	Violet
Tissue reaction	Minimal
Absorption	The hydrolytic action by which the material is broken down results in total absorption in approximately 180 to 210 days. Approximately 50% of tensile strength remain after 42 days.
Presentation	Box with 12 sutures



*For quality
&
safety*



www.sutures.be

SURGICRYL[®]

MONOFILAMENT

POLYDIOXANONE

CHARACTERISTICS

- High and long-lasting tensile strength
- Synthetic
- Breaks down by hydrolysis
- Monofilament with microscopic uniform diameter
- Pliability and softness
- Hermetically sealed packing

BENEFITS

- Outstanding strength retention for extended wound support
- Minimal tissue reaction
- Predictable absorption rate
- Soft passage through the tissues without sawing, tissue drag and trauma
- Absence of capillarity (no support of bacterial growth)
- Excellent knot tie down
- Good handling properties
- Easy and secure knotting
- Guaranteed seal and product sterility

INDICATIONS General soft tissue closing especially in


- General surgery
- Gastrointestinal surgery
- Ophthalmic surgery
- Orthopaedics
- Abdominal wall closure
- Gynaecology
- Plastic surgery
- Urology

CONTRAINDICATIONS Cardiovascular surgery
Neurosurgery
Prosthetic devices


Being absorbable not to use when wound support for more than 6 weeks is required.



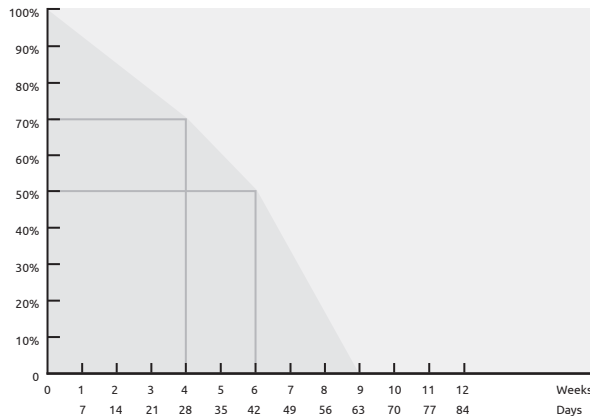


Caution in patients with delayed wound healing (e.g. compromised immune system).



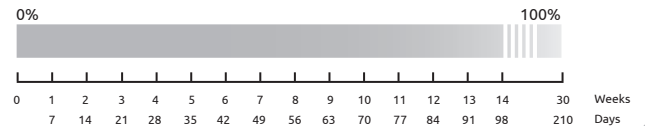






TENSILE STRENGTH



ABSORPTION



Round bodied taper point 		Metric	USP	Length		Order number
1/2 Circle 	HR 13 mm 	0.7	6/0	75 cm	<i>violet</i>	13070113
		1	5/0	75 cm	<i>violet</i>	13100113
	HR 17 mm 	0.7	6/0	75 cm	<i>violet</i>	13070117
		1	5/0	75 cm	<i>violet</i>	13100117
		1.5	4/0	75 cm	<i>violet</i>	13150117
		2	3/0	75 cm	<i>violet</i>	13200117
	HR 22 mm 	1.5	4/0	75 cm	<i>violet</i>	13150122
		2	3/0	75 cm	<i>violet</i>	13200122

Reverse cutting 		Metric	USP	Length		Order number
3/8 Circle 	DS 12 mm 	0.7	6/0	45 cm	<i>violet</i>	13071512
		1	5/0	45 cm	<i>violet</i>	13101512
		1.5	4/0	75 cm	<i>violet</i>	13151512
	DS 16 mm 	1	5/0	75 cm	<i>violet</i>	13101516
		1.5	4/0	75 cm	<i>violet</i>	13151516
	DS 16 mm black 	1	5/0	75 cm	<i>violet</i>	13101516VB
	DS 19 mm 	1	5/0	45 cm	<i>violet</i>	13101519
		1	5/0	75 cm	<i>violet</i>	13101519-75
		1.5	4/0	75 cm	<i>violet</i>	13151519
		2	3/0	75 cm	<i>violet</i>	13201519
	DS 19 mm Easy Pass 	1	5/0	50 cm	<i>violet</i>	13101516P
	DS 19 mm black 	1.5	4/0	75 cm	<i>violet</i>	13151519VB
	DS 24 mm 	1.5	4/0	75 cm	<i>violet</i>	13151524
		2	3/0	75 cm	<i>violet</i>	13201524
		3	2/0	75 cm	<i>violet</i>	13301524