

CATGUT

PLAIN



Type	Twisted multifilament with a monofilament appearance
Composition	Strands of purified collagen taken from the serosal layer of selected bovines. Origin from BSE-free classified countries.
Coating	None
Colour	Ivory
Tissue reaction	Moderate
Absorption	By phagocytosis. Complete mass absorption in approximately 63 days. Approximately 50% of tensile strength remain after 7 days. When used in infected tissues or in tissues with increased levels of proteolytic enzymes (stomach, cervix, vagina) it is more quickly absorbed.
Presentation	Box with 12 sutures



*For quality
&
safety*

smi

www.sutures.be

CATGUT

PLAIN

CHARACTERISTICS

- High tensile strength
- Absorption by phagocytosis
- Twisted structure

- Uniform diameter by centerless polishing (monofilament appearance)

- Hermetically sealed packing

BENEFITS

- Excellent tissue support
- Uniform absorption
- Excellent handling properties
- High flexibility
- High tensile strength
- Excellent knot security

- Soft passage through the tissues without sawing, tissue drag and trauma
- Excellent knot tie down

- Guaranteed seal and product sterility

INDICATIONS General soft tissue closing and/or ligation especially in

- General surgery
- Gastrointestinal surgery
- Gynaecology
- Obstetrics
- Urology
- Ophthalmic surgery

CONTRAINDICATIONS Cardiovascular surgery
Neurosurgery


Being absorbable not to use when an extended wound support is required.

Caution in patients with delayed wound healing (e.g. compromised immune system).


In the region of infected wounds

Known sensitivities or allergies to collagen

TENSILE STRENGTH

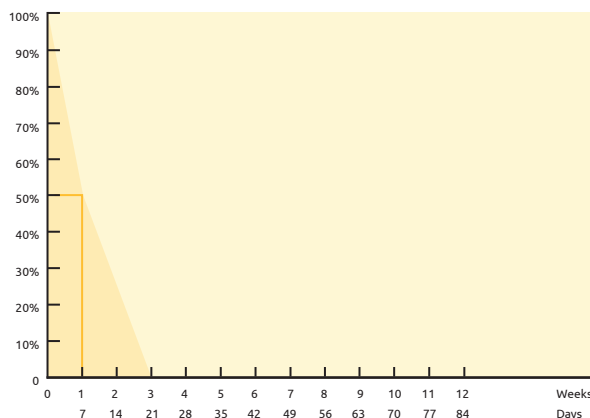

ABSORPTION

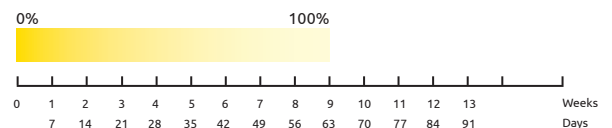

CATGUT

PLAIN

Without needle



	Metric	USP	Length		Order number
	2	4/0	150 cm	<i>ivory</i>	1200015
	3	3/0	150 cm	<i>ivory</i>	1300015
	3.5	2/0	150 cm	<i>ivory</i>	1350015
	4	0	150 cm	<i>ivory</i>	1400015
	5	1	150 cm	<i>ivory</i>	1500015
	6	2	150 cm	<i>ivory</i>	1600015
	7	3	150 cm	<i>ivory</i>	1700015

Round bodied taper point










	Metric	USP	Length		Order number
HR 16 mm	1.5	5/0	75 cm	<i>ivory</i>	1150116
	2	4/0	75 cm	<i>ivory</i>	1200116
HR 20 mm	2	4/0	75 cm	<i>ivory</i>	1200120
	3	3/0	75 cm	<i>ivory</i>	1300120
HR 26 mm	3.5	2/0	75 cm	<i>ivory</i>	1350120
	2	4/0	75 cm	<i>ivory</i>	1200126
HR 30 mm	3	3/0	75 cm	<i>ivory</i>	1300126
	3.5	2/0	75 cm	<i>ivory</i>	1350126
	4	0	75 cm	<i>ivory</i>	1400126
HR 37 mm	3	3/0	75 cm	<i>ivory</i>	1300137
	3.5	2/0	75 cm	<i>ivory</i>	1350137
	4	0	75 cm	<i>ivory</i>	1400137
	5	1	75 cm	<i>ivory</i>	1500137
HR 40 mm	5	1	90 cm	<i>ivory</i>	1500140-90
HR 45 mm	4	0	75 cm	<i>ivory</i>	1400145
	5	1	75 cm	<i>ivory</i>	1500145
	6	2	75 cm	<i>ivory</i>	1600145
HR 50 mm	3.5	2/0	75 cm	<i>ivory</i>	1350150

1/2 Circle



CATGUT

PLAIN

Reverse cutting 		Metric	USP	Length	Order number
3/8 Circle 	DS 12 mm 	1.5	5/0	45 cm	<i>ivory</i> 1151512
		1	6/0	45 cm	<i>ivory</i> 1101516
	DS 16 mm 	1.5	5/0	75 cm	<i>ivory</i> 1151516
		2	4/0	75 cm	<i>ivory</i> 1201516
		3	3/0	75 cm	<i>ivory</i> 1301516
	DS 19 mm 	2	4/0	75 cm	<i>ivory</i> 1201519
		3	3/0	75 cm	<i>ivory</i> 1301519
	DS 24 mm 	3.5	2/0	75 cm	<i>ivory</i> 1351519
		2	4/0	75 cm	<i>ivory</i> 1201524
		3	3/0	75 cm	<i>ivory</i> 1301524
	DS 24 mm	3.5	2/0	75 cm	<i>ivory</i> 1351524
		4	0	75 cm	<i>ivory</i> 1401524
DS 37 mm 	3.5	2/0	75 cm	<i>ivory</i> 1351537	

IN CASSETTES

Metric	USP	Length	Order number
2	4/0	100 m	1020
3	3/0	100 m	1030
3.5	2/0	100 m	1035
4	0	100 m	1040
5	1	75 m	1050
6	2	50 m	1060
7	3	50 m	1070
8	4	30 m	1080
9	5	30 m	1090
10	6	30 m	1100

